

PIPETBOY *comfort*



Accurate pipetting at your fingertips

PIPETBOY *comfort*

Flexible, lightweight and always ready to use

Charger and stand – all in one

Either with or without a pipette, on the table or on the wall – the utmost flexibility and ease of use for your laboratory. Your **PIPETBOY** *comfort* is on stand-by and fully charged, whenever you need it.



Light as a feather

A lightweight 190 grams, the strain is taken out of even prolonged pipetting work. Hands, arms and shoulders remain fatigue-free, enhancing your comfort.



Ergonomically designed

The **PIPETBOY** *comfort* sits snugly in the hand naturally matching its shape. The rounded handle also fits perfectly in small hands. Pipetting has never been so comfortable.



Flexible and compatible

The silicone Pipette Mount guarantees a sure grip on practically all available plastic or glass pipettes in the 1 – 100 ml volume range.



Maximum Safety

The easy-to-change sterile filter provides your **PIPETBOY** *comfort* with reliable protection from liquid contamination and damage. And if you do happen to aspirate critical liquids with the pipette, the **PIPETBOY** *comfort* silicone Pipette Mount, the Pipettor Nose Piece and the filter rubber are all autoclavable while the sterile filter can be replaced. So easy and so safe.



royal



classic



volcano



pacific



sahara




Battery indicator

An illuminated yellow LED shows when the battery needs to be charged. While the battery is being charged, the LED flashes to show that the contact has been correctly established. You know the status of your rechargeable battery at all times. When the **PIPETBOY comfort** is connected directly to the mains adapter, you can still pipette while it is charging. The **PIPETBOY comfort** is always on standby!



From rapid pipetting to sensitive precision dosing

You can control the speed of pipetting continuously via the knurled screw. Precise dosing is also possible by applying different amounts of pressure to the triggers.

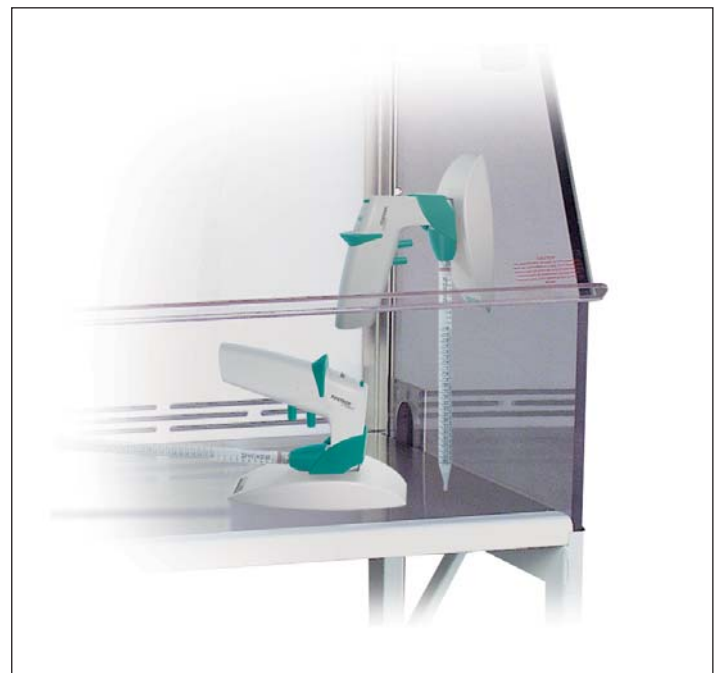

Held perfectly in all positions

Whether on the wall or tabletop, inside or outside a hood, the **PIPETBOY comfort** is held perfectly – with or without a pipette.

Wall mounting:

Magnet, adhesive strip or screw – the charger can be easily and securely mounted on any wall. The aerodynamic shape of the charger ensures minimal air turbulence in the laminar flow and reduces the risk of contamination.

Tabletop mounting:

Simply place the charger stand on top of your table – 4 rubber feet ensure stability.



Powerful

2 NiMH rechargeable batteries deliver optimum performance ensuring maximum power and durability with minimal weight. Fully charged batteries provide 8 hours of continuous operation. The **PIPETBOY comfort** can either be connected directly to the mains adapter or recharged via the charging console.

PIPETBOY comfort

Ordering information


| Description | Unit | Part No. |
|---|------|----------|
| PIPETBOY comfort classic , incl. charging station, adapter, sterile filter, charger mounting kit and two batteries | 1 | 156 000 |
| PIPETBOY comfort volcano , as Art. No. 156 000 | 1 | 156 001 |
| PIPETBOY comfort pacific , as Art. No. 156 000 | 1 | 156 002 |
| PIPETBOY comfort sahara , as Art. No. 156 000 | 1 | 156 003 |
| PIPETBOY comfort royal , as Art. No. 156 000 | 1 | 156 004 |

Ordering information for single and spare parts

| No. | Description | Unit | Part No. |
|-----|---|-------------------------|----------|
| 1 | Filter rubber | 1 | 153 225 |
| 2 | Filter sterile, 0.45 µm | 1 | 153 015 |
| | Filter unsterile, 0.45 µm | 1 | 153 016 |
| | Filter sterile, 0.2 µm | 1 | 153 608 |
| | Filter unsterile, 0.2 µm | 1 | 153 607 |
| 3 | Pipette Mount, silicone | 1 | 151 020 |
| 4 | Pipettor Nose Piece (classic) | 1 | 156 100 |
| | Pipettor Nose Piece (volcano) | 1 | 156 101 |
| | Pipettor Nose Piece (pacific) | 1 | 156 102 |
| | Pipettor Nose Piece (sahara) | 1 | 156 103 |
| | Pipettor Nose Piece (royal) | 1 | 156 104 |
| 5 | Battery-Set, 2 Pieces (1.2 V AAA NiMH/1000 mAh) | 2 | 156 151 |
| 6 | Battery cover Set (classic) | 1 | 156 110 |
| | Battery cover Set (volcano) | 1 | 156 111 |
| | Battery cover Set (pacific) | 1 | 156 112 |
| | Battery cover Set (sahara) | 1 | 156 113 |
| | Battery cover Set (royal) | 1 | 156 114 |
| 7 | Charging station (classic) | 1 | 156 050 |
| | Charging station (volcano) | 1 | 156 051 |
| | Charging station (pacific) | 1 | 156 052 |
| | Charging station (sahara) | 1 | 156 053 |
| | Charging station (royal) | 1 | 156 054 |
| 8 | Charger mounting kit (tape, magnet, base pads, screw) | 1 | 156 150 |
| 9 | Adapter: | | |
| | US version | 1 | 156 630 |
| | (100 – 240 VAC, 50/60 Hz) | | |
| | EU version | 1 | 156 631 |
| | UK version | 1 | 156 632 |
| | AU version | 1 </td <td>156 633</td> | 156 633 |
| | JP version | 1 | 156 634 |

Technical Data

| | | |
|-----------------------|---|-------------------------|
| Dimensions (LxHxW): | PIPETBOY comfort: | 130 x 145 x 35 mm |
| | Charging station: | 175 x 75 x 80 mm |
| | Mains adapter: | 75 x 45 x 45 mm |
| Weight: | PIPETBOY comfort: | 190 g |
| | Charging station: | 240 g |
| | Mains adapter: | 190 – 300 g |
| Batteries: | rechargeable, 2 x 1.2 V AAA NiMH/1000 mAh | |
| Mains adapter output: | 9 VDC/500 mA | |
| Operating condition: | Operation: | 5 – 40°C, max. 80% RH |
| | Storage: | -10 – 50°C, max. 95% RH |

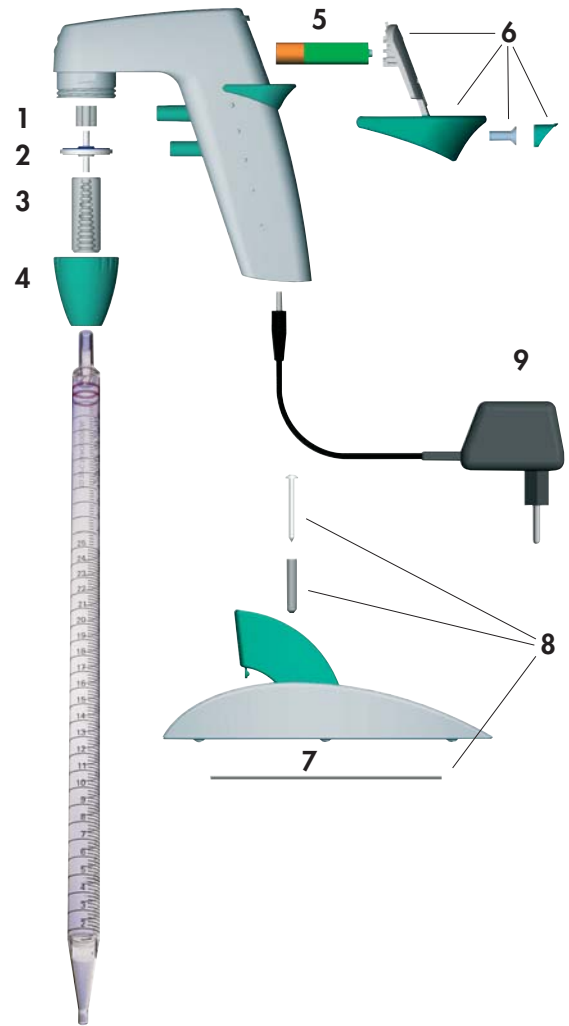

INTEGRA Biosciences AG

CH-7000 Chur, Switzerland
 Phone: +41 81 286 95 30
 Fax: +41 81 286 95 33
 info@integra-biosciences.com
 www.integra-biosciences.com